# Hypervolt

Operating Instructions | 取扱説明書

# EN

## PLEASE READ PRIOR TO USE

The Hypervolt is a handheld percussion device that helps your warm up and feel better even quicker before and after physical activity.

## OPERATING

Switch your Hypervolt to the ON position at the base of the handle. Once your device is powered ON a blue light will illuminate, plus (up to) four green lights for battery level. Next, activate the pressure sensor by pressing the speed setting button once. Wait five seconds as the white lights sequence and turn off. Your Hypervolt is now ready for use.

Change speeds (1-3) by pressing the speed setting button once per speed. The blue lights above the speed setting button will indicate the current speed.


To change head attachments insert the desired head directly into the shaft / opening while pushing firmly. For the Fork attachment align the "groove" on the fork with the recess on the shaft / opening and push firmly.

## CLEANING AND STORAGE

Make sure the power is OFF and the battery charger is not attached. Use a damp, clean cloth and gently wipe your Hypervolt. Store in a clean, cool place when not in use.

## SPECIFICATIONS

1. Removable battery release
2. Handle / rechargeable lithium-ion battery
3. Speed indicators (3)
4. Speed setting button
5. Pressure sensor indicators (3)
6. Battery level indicators (4 / LED band)
7. ON / OFF power switch
8. Interchangeable head attachments (5)
9. Bluetooth® connection indicator

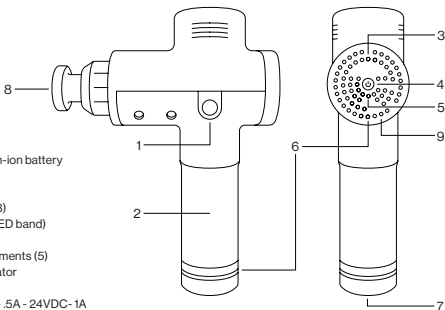AC charger: 100-240V 50/60Hz - .5A - 24VDC - 1A

Frequency: level 1 - 1800/30 Hz, level 2 - 2400/40 Hz, level 3 - 3180/53 Hz

Weight: 2.5 lbs/1.1 kg

**DANGER: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, OR PROPERTY DAMAGE, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING WARNINGS, CAUTIONS AND SAFETY INSTRUCTIONS:**

## **IMPORTANT SAFETY INSTRUCTIONS**

The Hypervolt (the term "Hypervolt" shall refer to both the Hypervolt and Hypervolt Plus models) is for adult use only and should always be used in accordance with the following warnings and instructions.

When using the Hypervolt use only on dry, clean surfaces of the body by lightly pressing and moving across the body for approximately sixty(60) seconds per region. The Hypervolt should only be used on soft tissue areas of the body and should not be used on any hard or bony areas of the body, including the head. Please immediately stop use of the Hypervolt upon experiencing any pain or discomfort, apart from light muscle soreness. Bruising can occur from use of the Hypervolt, regardless of the pressure setting, and if bruising does occur you should not use the Hypervolt on the bruised area until the bruise is fully healed. Use of the Hypervolt should also be avoided on all bruises, contusions, rashes or irritated or injured areas of the skin until they are fully healed. When using the Hypervolt be sure to keep fingers, toes, hair, and other body parts away from the back of the applicator to avoid pinching or hair entanglement.

The Hypervolt should be kept dry, should never be immersed in water and should be kept away from liquids and heat sources at all times. The ventilation holes on the Hypervolt should be kept free of dust, debris and/or any other substance that could penetrate the holes, including but not limited to liquids, dirt, sand, or food. You should never operate the Hypervolt for more than one (1) continuous hour at a time and after each use please allow the Hypervolt to rest for a minimum of thirty (30) minutes before using it again.

**THE HYPERVOLT SHOULD ONLY BE CHARGED WITH THE TWO PRONGED, DOUBLE INSULATED, 26VDC HYPERICE CHARGER (MODEL NUMBER SK03T1-2600100W2 OR MX24D2-2601000) THAT WAS SUPPLIED WITH THE UNIT. DO NOT CHARGE THE HYPERVOLT OVERNIGHT OR LEAVE THE HYPERVOLT UNATTENDED WHEN CHARGING OR IN USE.**

After charging the Hypervolt, the charger should be unplugged from the Hypervolt's battery handle prior to use. Prior to disposing of a Hypervolt unit the battery must be removed. The battery must be complied with safety at an appropriate e-waste disposal or recycling facility.

The Hypervolt should never be used for any purpose other than its intended use, should not be dropped, and should never be tampered with or altered in any way.

**PLEASE DO NOT USE THE HYPERVOLT, OR ANY PERCUSSION OR VIBRATION DEVICE WITHOUT FIRST OBTAINING APPROVAL FROM YOUR PHYSICIAN IF ANY OF THE FOLLOWING APPLY:**

Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pace-makers, recent surgery or injury, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUD's, metal pins or plates or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any percussion or vibration device. These contra indications do not mean that you are not able to use a percussion or vibration device but we advise that before doing so you consult your physician. Ongoing research is being done on the effects of percussion massage for specific medical disorders that may result in the shortening of the list of contra indications as shown above.

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Copyright 2020 © Hyperice, Inc. All rights reserved.